

Supercath™ AZ

Safety Dialysis Catheter

Ultimate Safety Technology

For Dialysis Patients and Health Care Workers




Latest 'Safety' design with fully-encased needle, **Supercath™ AZ** helps to prevent needlestick injury and to minimize risk of blood exposure, provides passive safety technology to dialysis patients and health care workers.



Supercath™ AZ

Safety Dialysis Catheter

- Features and Benefits -

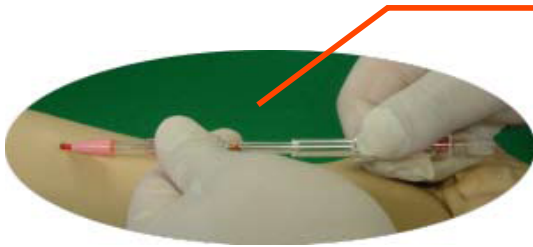

1) A Unique, Easy Handling

Telescopic Needle Retraction System

Eliminates accidental needlestick for both health care workers and patients on removal of introducer needle, and prevents blood exposure with fully-encased needle



2) FEP HD Catheter

- Minimizes risks of internal AV fistulae damage
- Extremely comfortable compared to metal AVF needles
- Highly biocompatible



3) Male Cap with Integrated Stop Valve

No blood exposure after removal of introducer needle

4) Ergonomically Tapered Catheter

Catheter tapered up to larger hub diameter enhances flow



5) Lancet Point Needle Tip

Minimizes vessel trauma with smooth cannulation and lessens patients' discomfort

6) Single Piece Body Molding Process

Reduces hemolysis with maximized flow

Product No.	Catheter (mm)			Inner Needle Gauge	No. of Side Holes	Color Code	Packing
	Length	O.D.	I.D.				
SP302-16(25)	25	2.10 (14G)	1.60 (16G)	16G	4	Green	100pcs/Shelf Box 1000pcs/Shipper Box (10 Shelf Boxes)
SP302-16(33)	33						
SP302-17(25)	25	1.93 (15G)	1.45 (17G)	17G	4	Orange	
SP302-17(33)	33						
SP302-18(25)	25	1.73 (16G)	1.25 (18G)	18G	4	Pink	
SP302-18(33)	33						

PRODUCT NAME: SUPERCATH SAFETY CLAMPCATH



Distributed exclusively in Canada by



199 Saint David Street
Lindsay, Ontario, Canada
K9V 5K7

Tel: 705.328.2518
Fax: 705.328.9747
Toll-Free: 800.387.9757
mail@cardiomed.com

Site: www.cardiomed.com

Manufacturer: Togo Medikit Co., Ltd.

HDBA080001-A60F2H-A