

Peritoneal dialysis catheters

Performance and comfort

Available in a variety of styles CardioMed's Peritoneal dialysis catheters offer consistent performance and patient comfort.

CardioMed's Peritoneal dialysis catheters offer the following features:

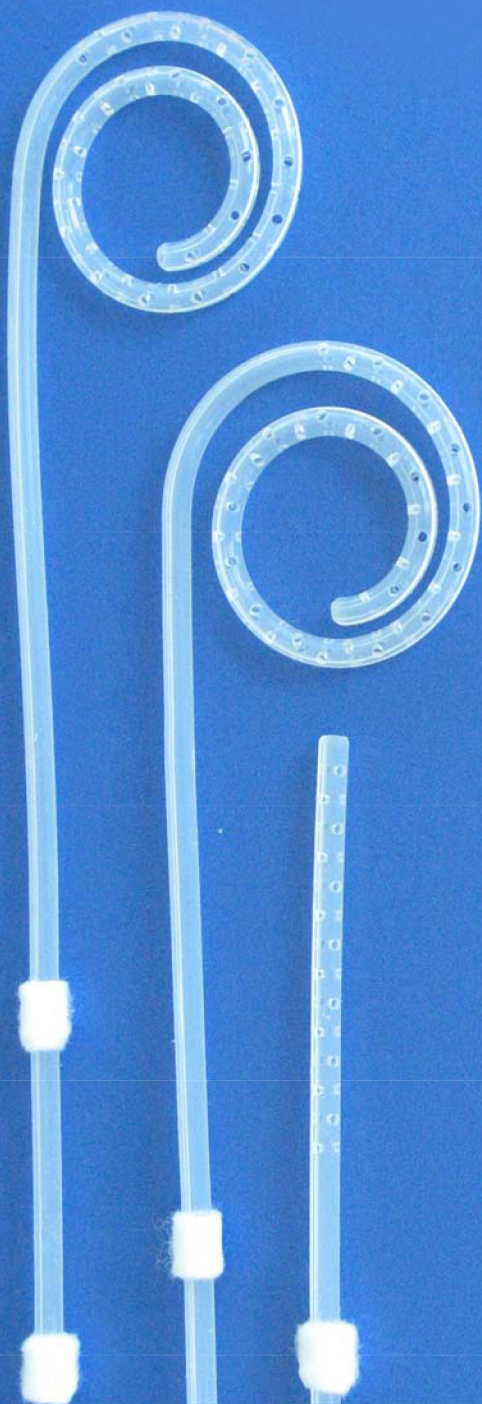
Larger catheter ID lumen and larger drainage holes allow for faster inflow and outflow rates

Radiopaque stripe enables easy visualization of the catheters position under X-Ray

Precise positioning of holes ensures that a solid continuous radiopaque stripe is visualized

Dacron cuff(s) facilitate aggressive tissue ingrowths

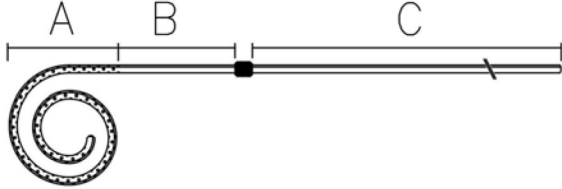
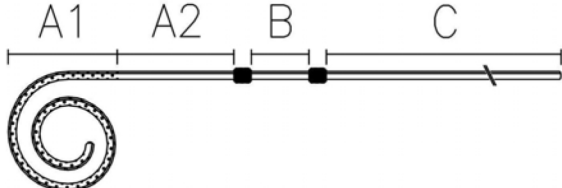
Beveled tip reduces trauma on abdominal tissues



CardioMed
Life-sustaining solutions

Peritoneal dialysis catheters

ORDERING INFORMATION

SINGLE CUFF STRAIGHT CATHETERS		LENGTHS (cm)				PROFILE	
ITEM	DESCRIPTION	TOTAL	A	B	C		
S4-120-SSP	Straight peritoneal dialysis catheter	20	5.1	5.9	8.0		
S4-124-SSP	Straight peritoneal dialysis catheter	24	7.0	7.1	8.9		
S4-132-SSP	Straight peritoneal dialysis catheter	32	8.4	8.7	13.9		
S4-142-SSP*	Straight peritoneal dialysis catheter	42	10.3	12.4	18.3		
DOUBLE CUFF STRAIGHT CATHETERS		LENGTHS (cm)					
ITEM	DESCRIPTION	TOTAL	A1/A2	B	C		
S4-232-SSP	Straight peritoneal dialysis catheter	32	7.0/3.2	7.1	12.7		
S4-242-SSP	Straight peritoneal dialysis catheter	42	10.2/5.1	5.0	19.7		
SINGLE CUFF COILED CATHETERS		LENGTHS (cm)					
ITEM	DESCRIPTION	TOTAL	A	B	C		
S4-155-CSP*	Coiled peritoneal dialysis catheter	64	19.0	16.5	27.5		
DOUBLE CUFF COILED CATHETERS		LENGTHS (cm)					
ITEM	DESCRIPTION	TOTAL	A1/A2	B	C		
S4-255-CSP*	Coiled peritoneal dialysis catheter	64	19.5/16.5	5.0	21.0		

Except where otherwise noted, all products are latex-free

* To order catheters with an IR insertion kit, add IR to the end of the product code

Unless otherwise indicated, products listed are ETO sterilized. Sets are packaged 10 per box. IR kits are 5 per box.

A Peritoneal catheter set includes:

- 1 catheter
- 2 sealing caps
- 1 female luer adaptor
- 1 high viscosity gel pack

A Peritoneal catheter, IR insertion kit includes:

- 1 catheter
- .038"x70cm J guidewire
- 18Gx2 ¾" seldinger needle
- 16F peel-away sheath and dilator
- 2 sealing caps
- 1 female luer adaptor
- 1 high viscosity gel pack



CardioMed Supplies Inc.

199 Saint David Street,
Lindsay, ON K9V 5K7 Canada
Tel: 705.328.2518 Toll-free: 800.387.9757
Fax: 705.328.9747
Web: www.cardiomed.com Email: mail@cardiomed.com