


INTRAVASCULAR THERAPIES

PRESSURE INJECTABLE
IMPLANTABLE INFUSION PORTS

polysite

Comfort your patient
deserves. Strength you
demand.

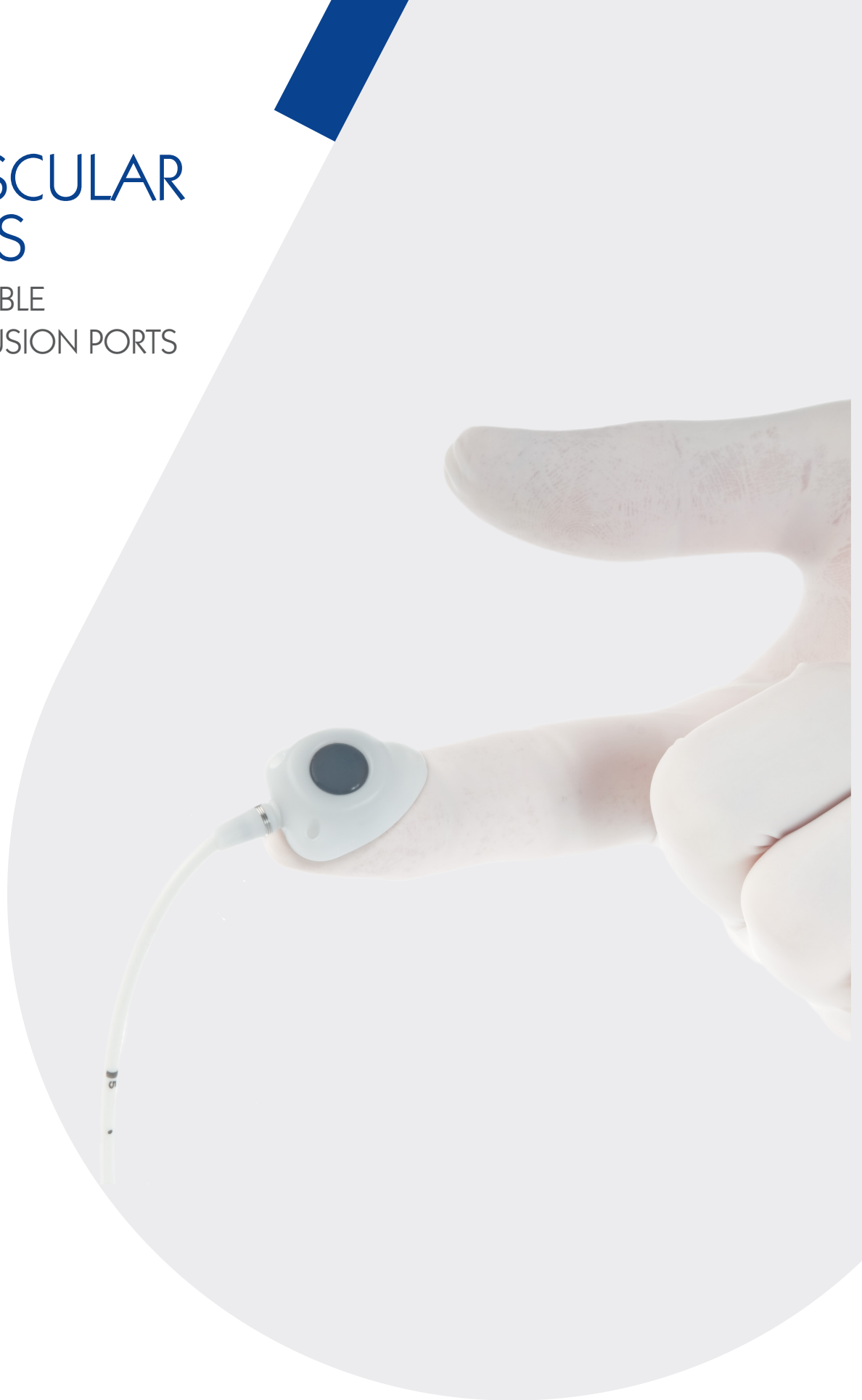


Value Life

polysite hybrid ports

Comfort your patient deserves.

Port combines lightweight plastic body for comfort and a strong titanium reservoir for durability

- Pressure Injectable up to 5 mL/sec at 300 psi*
- Silicone-filled suture holes available for easy explantation
- High flow rates

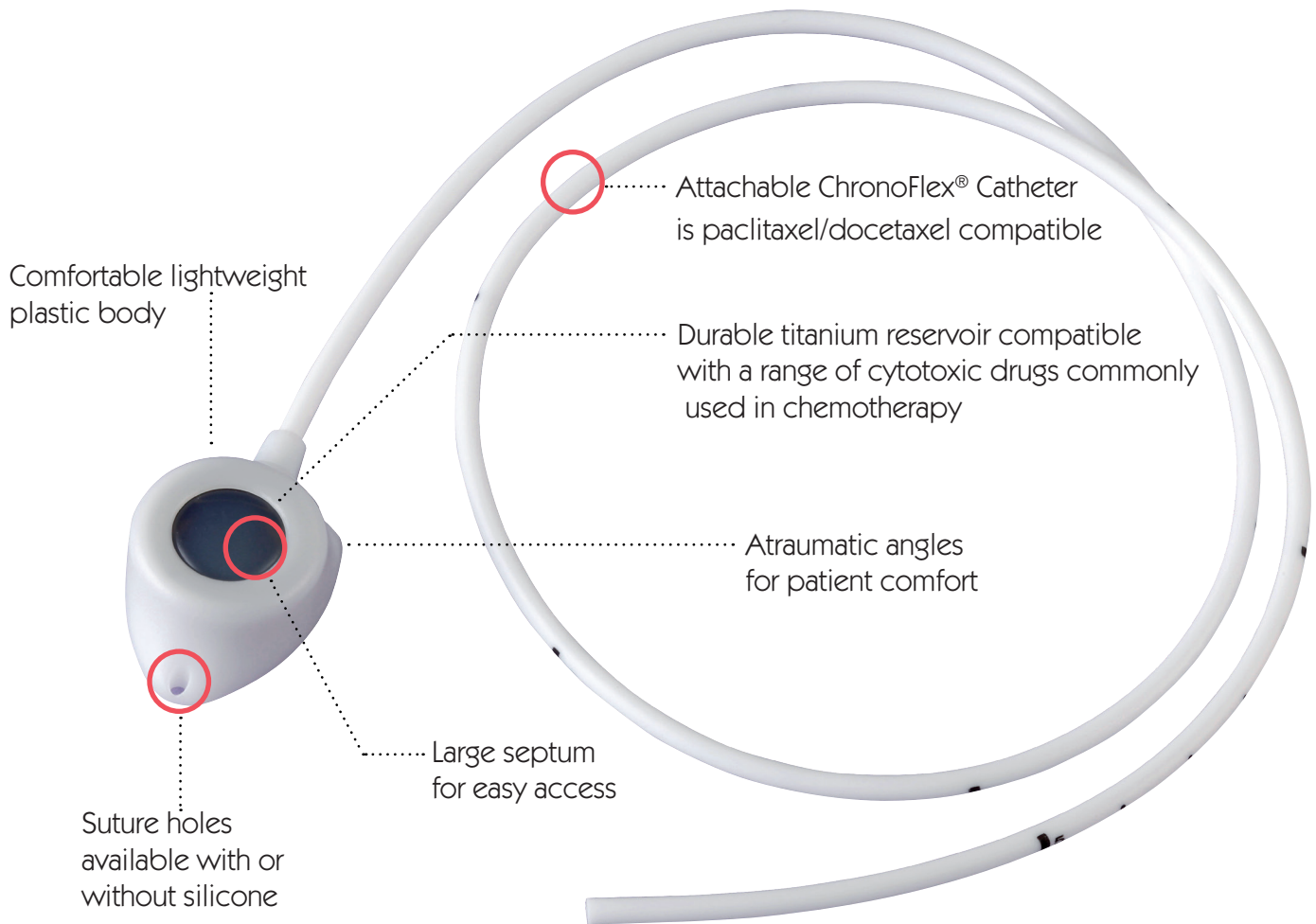


Available in 6F Micro, 7F Mini, 8F Standard sizes

polysite hybrid ports

Strength you demand.



MR Conditional

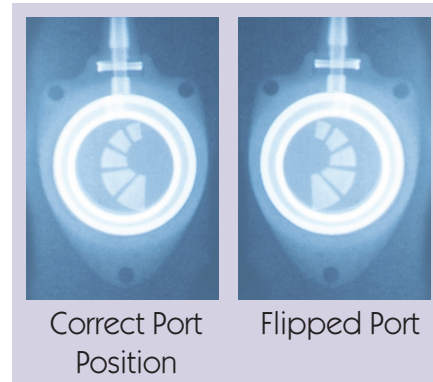

Correct Port Position

Flipped Port Position


Pressure Injectable*

- Easy to identify radiopaque power symbol
- Pressure Injectable up to 5 mL/sec at 300 psi*
- MR conditional at 1.5 and 3 Tesla

Features

- Catheter attaches easily to stem with a radiopaque locking collar
- Locking collar extends past port stem to reduce risk of stress on the catheter
- Internal diameter of port stem increases with catheter French size to enhance flow rate

Safety for you and your patients

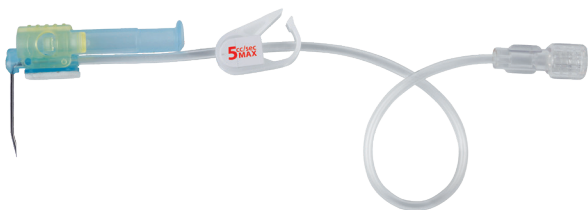

Positive Pressure Safety Huber Needle

- Patented needle removal system maintains positive flush pressure during one-handed needle removal
- 5 mL/sec infusion at 300 psi



Peelable Valved Sheath

- Reduce risk of air embolism or blood loss
- Smooth dilator to sheath transition

*Pressure injectable models only

Refer to package insert for warnings, precautions, contraindications, potential adverse effects, MR Conditional information and Instructions For Use. RX Only

Ordering Information

POLYSITE® Ports Micro 2000 Series

Product Code	Description	Silicone Filled Suture Holes	With Saftey Sheath	Pressure Injectable (5 mL/sec at 300 psi)	French Size	Catheter ID (mm)	Flow Rate (mL/H)*	Units/Box
2016SPI	POLYSITE® Pressure Injectable Implantable Infusion Port with 6F 45cm Attachable ChronoFlex® Catheter	•	•	•	6	1.30	1327	1
2016PI	POLYSITE® Pressure Injectable Implantable Infusion Port with 6F 45cm Attachable ChronoFlex® Catheter		•	•	6	1.30	1327	1

*Flow Rate (mL/h) using 19Ga Huber Needle

POLYSITE® Ports Mini 3000 Series

Product Code	Description	Silicone Filled Suture Holes	With Saftey Sheath	Pressure Injectable (5 mL/sec at 300 psi)	French Size	Catheter ID (mm)	Flow Rate (mL/H)*	Units/Box
3017SPI	POLYSITE® Pressure Injectable Implantable Infusion Port with 7F 60cm Attachable ChronoFlex® Catheter	•	•	•	7	1.45	1581	1
3017PI	POLYSITE® Pressure Injectable Implantable Infusion Port with 7F 60cm Attachable ChronoFlex® Catheter		•	•	7	1.45	1581	1

*Flow Rate (mL/h) using 19Ga Huber Needle

POLYSITE® Ports Standard 4000 Series

Product Code	Description	Silicone Filled Suture Holes	With Saftey Sheath	Pressure Injectable (5 mL/sec at 300 psi)	French Size	Catheter ID (mm)	Flow Rate (mL/H)*	Units/Box
4018SPI	POLYSITE® Pressure Injectable Implantable Infusion Port with 8F 60cm Attachable ChronoFlex® Catheter	•	•	•	8	1.70	2042	1
4018PI	POLYSITE® Pressure Injectable Implantable Infusion Port with 8F 60cm Attachable ChronoFlex® Catheter		•	•	8	1.70	2042	1

*Flow Rate (mL/h) using 19Ga Huber Needle

Pressure Injectable Kit Contents Include:

- One Pressure Injectable Low Profile Hybrid Port with Silicone Septum
- One Radiopaque ChronoFlex® Polyurethane Catheter
- One Straight Huber Needle 22Ga
- One Polyperf Safe™ Safety Huber Needle with Y-Site 20Ga x 25mm
- One Peelable Valved Introducer with Tissue Dilator
- Two Catheter Connection Rings with Radiopaque Marker Bands
- One 0.035cm x 60cm Marked Guidewire with "J" Tip on One and Straight Soft Tip on the other
- One Introducer Needle 18Ga x 70mm
- One Tunneler
- One 10mL Syringe
- One Vein Pick
- One Flushing Connector

FOR FURTHER INFORMATION, PLEASE CONTACT: customerservice-us@vygon.com

The specification shown in this leaflet are for information only and are not, under any circumstances, of a contractual nature.