


# GUIDE PLUS™ II

guide extension



**NIPRO**  
CANADA

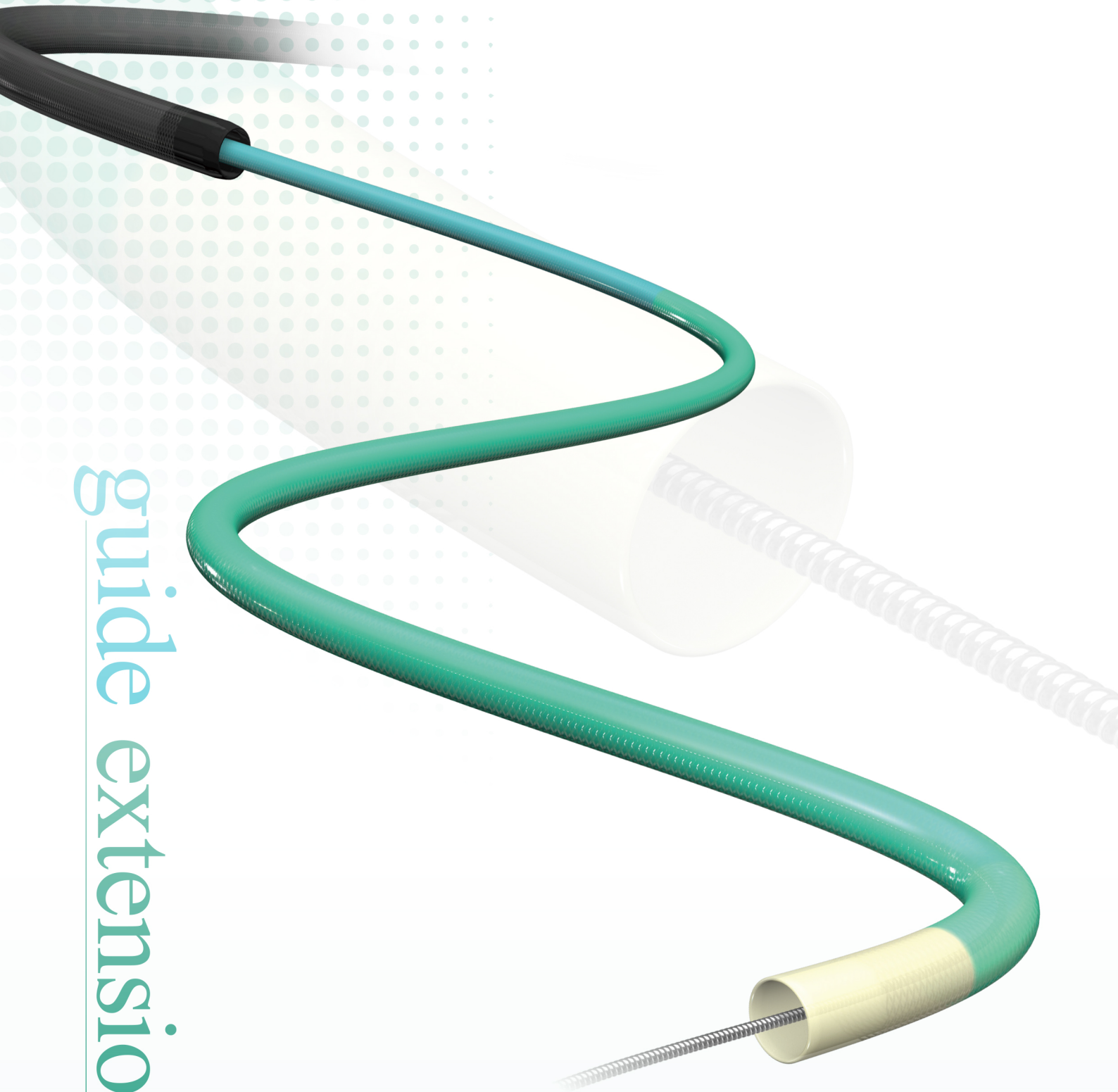




# Simplifying problematic cases by making delivery easier

---

A combination of flexible shaft and hydrophilic coating ensures excellent crossability in cases whereby direct stent delivery is problematic.

GUIDE  PLUS™ II

# FEATURES

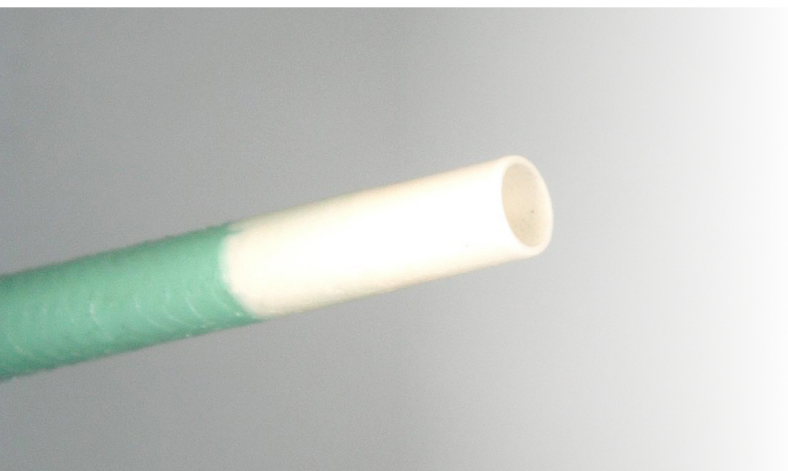

## Excellent flexible tip

The GUIDEPLUS II tip utilizes a flexible material that is rounded in shape to provide for safe delivery through the vasculature

## Excellent crossability

The braided tubing enhances pushability while the hydrophilic coating realizes excellent deliverability



## Smooth device insertion

A unique shaped port provides for easy device insertion


## Variation

2 types available:  
ST: Superior Trackability & EL: Expanded Lumen




# CATHETER DESIGN

## GUIDEPLUS II ST

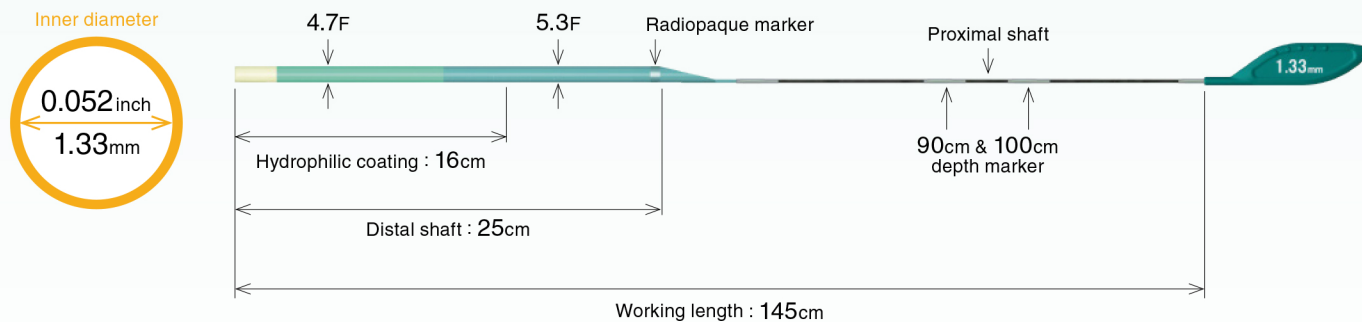

## GUIDEPLUS II EL

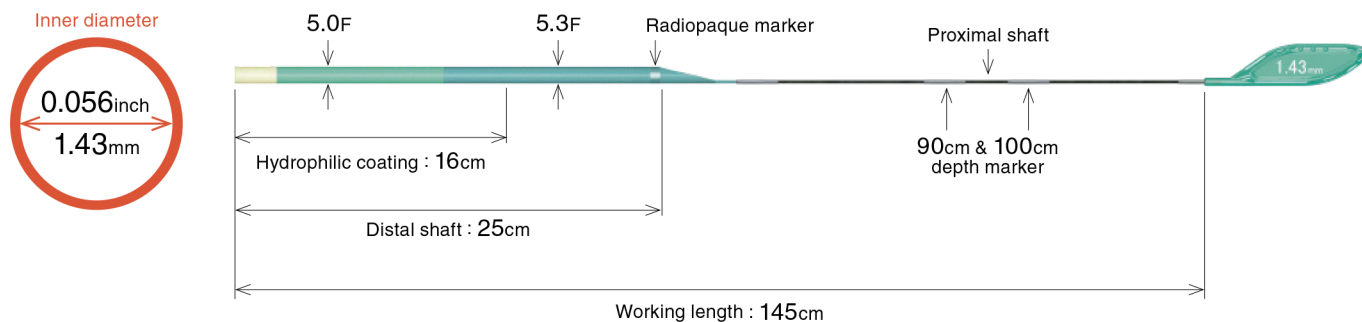

# ORDERING INFORMATION

Product name	Model no.	Max. external diameter	Distal shaft length	Compatible guiding catheter	Inner diameter
GUIDEPLUS II	GUIDEPLUS 2 ST INT	1.76mm / 5.3F	250mm	6F / $\geq 0.070$ inch	1.33mm / 0.052inch
	GUIDEPLUS 2 EL INT				1.43mm / 0.056inch



<https://www.nipro.ca>

**NIPRO CANADA**  
199 Saint David Street,  
Lindsay, Ontario, Canada  
T: 705-328-2518  
F: 705-328-9747  
E: service@nipro.ca

- ▶ Please consult your sales representative regarding stock availability.
- ▶ GUIDEPLUS is a trademark of NIPRO Corporation.
- ▶ Please refer to the Instructions For Use for further details regarding product use and specifications.

[Manufacturing facility]  
**NIPRO CORPORATION**  
3-26, Senriokashinmachi, Settsu, Osaka, 566-8510, Japan

GP11-2303-E9